The medical ProfiLine range



LKPv/LGPv: Range

Liebherr medical range 2007

The medical ProfiLine range LKPv/LGPv: Range



LKPv 6520

LKPv 6522 LKPv 1422

LKPv 1420

LGPv 6520 LGPv 1420

Liebherr laboratory refrigerators & freezers

Common features:

- High quality white enamelled exterior 304 grade stainless steel interior
- Powerful forced-air refrigeration system with rear air-duct & hot-gas defrost for minimum temperature variation
- Optical & acoustic hi-lo temperature and door-open alarms
- Volt-free contact for external alarm
- RS 485 interface: RS 232 adapter optional extra
- Optional infrared data transfer
- Integrated data logger & battery pack: saves temperature data for 72 hrs in case of power failure
- 83 mm insulation for maximum temperature stability
- High quality electronic controller with text display. Beside the pre-set temperature programs all temperature and alarm parameters can be adjusted to suit your individual application
- Control panel can be protected using a password to avoid unauthorised adjustments to the settings
- Castors fitted as standard
- Door lock, self-closing doors, fan cut-out switch





The medical ProfiLine range LKPv 6520 / LKPv 6522



LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

Single-door laboratory refrigerators

LKPv 6520 – solid door LKPv 6522 – glass door

Temperature range:

LKPv 6520, solid door: -2° C / +15 $^{\circ}$ C LKPv 6522, glass door: 0° C / +15 $^{\circ}$ C

- 4 coated wire shelves on adjustable trayslides, 60 kg shelf load
- Dimensions: 215 x 70 x 83 cm (H x W x D)
- 3 pre-set storage settings: blood, pharmacy, chromatography

Climate class:

LKPv 6520, solid door: +43 °C LKPv 6522, glass door: +32°C





The medical ProfiLine range LKPv 1420 / LKPv 1422



LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

Double-door laboratory refrigerators

LKPv 1420 – solid doors LKPv 1422 – glass doors

Temperature range:

LKPv 1420, solid door: -2° C / +15°C LKPv 1422, glass doors: 0° C / +15°C

- 8 coated wire shelves on adjustable trayslides, 60 kg shelf load
- Dimensions: 215 x 143 x 83 cm (H x W x D)
- 3 pre-set storage settings: blood, pharmacy, chromatography

Climate class:

LKPv 1420, solid door: +43 °C LKPv 1422, glass door: +32°C



The medical ProfiLine range LGPv 6520



LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

Single-door laboratory freezer

LGPv 6520 - solid door

Temperature range: -10°C / -35°C

- 4 coated wire shelves on adjustable trayslides, 60 kg shelf load
- Dimensions: 215 x 70 x 83 cm (H x W x D)
- 2 pre-set storage settings: laboratory freezer, blood plasma

Climate class LGPv 6520: +43 °C



The medical ProfiLine range LGPv 1420



LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420 **Double-door laboratory freezer**

LGPv 1420 - solid doors

Temperature range: -10°C / -35°C

- 8 coated wire shelves on adjustable trayslides, 60 kg shelf load
- Dimensions: 215 x 143 x 83 cm (H x W x D)
- 2 pre-set storage settings: laboratory freezer, blood plasma

Climate class: LGPv 1420: +32 °C



The ProfiLine evolution


LIEBHERR

LKPv/LGPv: access ports

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

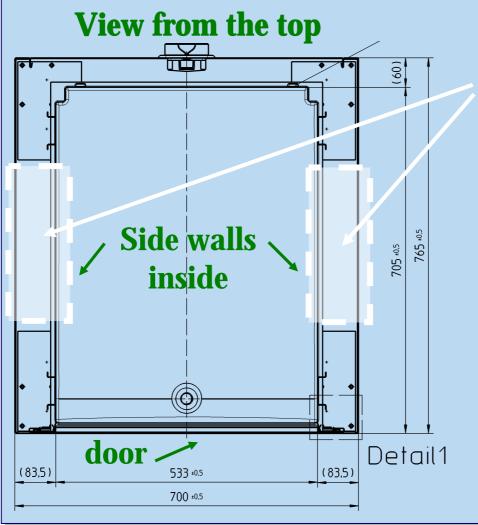

Ideal area for drilling through the sidewall.

Note: access ports are not available as a factory option.

Advantage of the flat inner liner

- easy to drill
- ✓ easy to seal

The ProfiLine evolution



LKPv/LGPv: Inner liner features

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

Coated wire shelves on adjustable trayslides:

- **✓** for general laboratory applications
- **✓** flexible use of storage volume
- **✓** non-tilting shelves

Smooth 304-grade s/s inner liner:

- ✓ hygienic, easy to clean & disinfect
- **✓** top standard
- ✓ acid resistant and solvents



The ProfiLine evolution



LKPv/LGPv: Inner liner features

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

robust s/s ladderracks

- ✓ small shelf spacing (27,5 mm) for ideal setting and adjusting of the shelves
- **✓** easily removable for cleaning

robust s/s "U"-trayslides

- ✓ shelves are easy to adjust
- **✓** shelves are non-tilting

Robust shelves

✓ max. load 60kg

Glass-door model with interior light

- ✓ Easy control of stored goods without opening the door
- **✓** Stored goods are illuminated for even better visibility

The ProfiLine evolution



Premium electronic controller

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420 • Temperature control: Premium electronic controller

controller with clear text display

• Temperature/text display: digital, illuminated

• Keyboard lock function: included

• data logger / records 3500 measurements

HACCP function:

• Temperature-recording: every 2 minutes

equals: ~ **116** hours / ~ **4,8** days

The ProfiLine evolution



Premium electronic controller

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420 • Temperature alarm,

high & low temp.: optical and acoustic

• Door open-alarm: optical and acoustic (4 min. delay)

• door open: displayed as text "door open"

• volt-free contact: included

• Interfaces: RS 485 & Infrared

• Battery to maintain functions of the controller

during power failure: saves data for 72 hours

The ProfiLine evolution



Premium electronic controller, pre-set programs

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

Pre-set programs Lab cooling

Pre-set Programs	Name	Relative Humidity	Set temperature
1	Blood	medium	+4°C
2	Pharmacy	low	+5°C
3	Chromatography	low	+4°C

The ProfiLine evolution



Premium electronic controller, pre-set programs

LGPv 6520 LGPv 1420

Pre-set programs Lab freezing

Pre-set Programs	Name	Relative Humidity	Set temperature
1	Lab	high	-20°C
2	Blood plasma	high	-30°C

The ProfiLine evolution



Premium electronic controller, accessories

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

GKPv 6580 GKPv 1480

GGPv 6580 GGPv 1480

standard interface adapter RS 485 to RS 232 with data logging software

Connection module from standard bus RS 485 to RS 232 to save data



- RS 232

Data transfer possible up to max. 10m

RS 485
Data transfer possible up to max. 1000 m

The ProfiLine evolution



electronic, Accessories

LKPv 6520 LKPv 1420

LKPv 6522 LKPv 1422

LGPv 6520 LGPv 1420

GKPv 6580 GKPv 1480

GGPv 6580 GGPv 1480 Accessory: infrared key with software

Transfer of data from controller to PC

With the IR key the data can be read-out and transferred to the PC



LKPv 6520

The ProfiLine evolution



LKPv/LGPv: Technical data

Technical data LKPv 6520 Lab forced- air refrigerator

Dimension exterior: Width: 700 / 533 mm interior: Depth: 830 / 650 mm

Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 2,0 A; 250 Watt; 220-230 V / 50 Hz

Energy consumption: 1,25 kWh/24h Temperature range: -2° to +15°C

Refrigerant: R 134a

Climate class: T (+18°C to +43°C)

Gross volume: 601 litre

shelves: 4 coated wire shelves

Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Insulation: 83 mm (door 60 mm)

Door: foamed Interior lighting: none

Castor: 2 fixed and 2 swivel with brake

The ProfiLine evolution



LKPv/LGPv: technical data

LKPv 6522

LKPv 6522 Lab forced- air refrigerator

Dimension exterior : Width: 700 / 533 mm interior: Depth: 830 / 650 mm Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 2,5 A; 300 Watt; 220-230 V / 50 Hz

Energy consumption: 2,7 kWh/24h Temperature range: 0° to +15°C

Refrigerant: R 134a

Climate class: $N (+16^{\circ}C \text{ to } +32^{\circ}C)$

Gross volume: 601 litre

shelves: 4 coated wire shelves

Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Isolation: 83 mm (door 60 mm)

Door: glass door

Interior lighting: Illuminated 15 Watt

Castor: 2 fixed and 2 swivel with brake

The ProfiLine evolution



LKPv/LGPv: technical data

LKPv 1420

LKPv 1420 Lab forced- air refrigerator

Dimension exterior : Width: 1430 / 1263 mm interior: Depth: 830 / 650 mm Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 2,5 A; 350 Watt; 220-230 V / 50 Hz

Energy consumption: 1,9 kWh/24h Temperature range: -2° to +15°C

Refrigerant: R 134a

Climate class: T (+18°C to +43°C)

Gross volume: 1427 litre

shelves: 8 coated wire shelves

Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Isolation: 83 mm (door 60 mm)

Door: foamed Interior lighting: none

Castor: 3 fixed and 2 swivel with brake

The ProfiLine evolution



LKPv/LGPv: Technical data

LKPv 1422

LKPv 1422 Lab forced- air refrigerator

Dimension exterior : Width: 1430 / 1263 mm interior: Depth: 830 / 650 mm Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 3,0 A; 400 Watt; 220-230 V / 50 Hz

Energy consumption: 5,0 kWh/24h Temperature range: 0° to +15°C

Refrigerant: R 134a

Climate class: $N (+16^{\circ}C \text{ to } +32^{\circ}C)$

Gross volume: 1427 litre

shelves: 8 coated wire shelves

Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Isolation: 83 mm (door 60 mm)

Door: glass door

Interior lighting: Illuminated 2 x 15 Watt

Castor: 3 fixed and 2 swivel with brake

The ProfiLine evolution



LKPv/LGPv: Technische Daten

LGPv 6520

LGPv 6520 Lab forced- air freezer

Dimension exterior: Width: 700 / 533 mm interior: Depth: 830 / 650 mm

Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 4,0 A; 600 Watt; 220-230 V / 50 Hz

Energy consumption: 5,5 kWh/24h Temperature range: -10° to -35°C

Refrigerant: R 404a

Climate class: T (+18°C to +43°C)

Gross volume: 601 litre s/s grids: 4 pieces Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Isolation: 83 mm (door 60 mm)

Door: foamed Interior lighting: none

Castor: 2 fixed and 2 swivel with brake

The ProfiLine evolution



LKPv/LGPv: Technische Daten

LGPv 1420

LGPv 1420 Lab forced- air freezer

Dimension exterior: Width: 1430 / 1263 mm interior: Depth: 830 / 650 mm

Height: 2150 / 1550 mm

Cabinet: steel sheet 0,7 mm powder coated

Inner liner: s/s (AISI 304) 0,5 mm

Connection rating: 4,5 A; 700 Watt; 220-230 V / 50 Hz

Energy consumption: 9,8 kWh/24h Temperature range: -10° to -35°C

Refrigerant: R 404a

Climate class: $N (+16^{\circ}C \text{ to } +32^{\circ}C)$

Gross volume: 1427 litre

shelves: 8 coated wire shelves

Lock: standard

Control: Electronic with clear text, HACCP function

Defrost: Hot-gas

Isolation: 83 mm (door 60 mm)

Door: foamed Interior lighting: none

Castor: 3 fixed and 2 swivel with brake